

# **Disposable Dental Air Water Syringe Tips**




AWT01,universal three-in-one Air Water Syringe Tips. The inner tube is made of stainless steel and the outer tube is made of white plastic.

Specification: 84X  $\varphi$  3.87mm Packing: 250 pcs/box,30 boxes/ctn Carton Size: 290X260X230mm Weight: 11.6 KGS

White \_\_\_



# **Disposable Dental Air Water Syringe Tips**

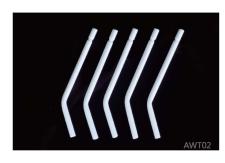


AWT02, white plastic material. Independent plastic tube structure without core. There are two modes with Syringe to use: water spray and water

Item: AWT02 Specification: 80X  $\varphi$  3.87mm Packing: 250 pcs/box,30 boxes/ctn Carton Size: 290X260X230mm Weight: 7 KGS



# **Disposable Dental Air Water Syringe Tips**

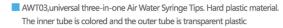

White



Item: AWT03 Specification: 92X \$\dagger\$ 3.87mm Packing: 250 pcs/box,40boxes/ctn Carton Size: 420X270X230mm Weight: 13.1 KGS



# **Disposable Dental Air Water Syringe Tips**



Item: AWT05 Specification: 76X \$\Phi\$ 3.87mm Packing: 250 pcs/box,40 boxes/ctn Carton Size: 340X270X230mm



# **DAM KIT**







### **Dental Dam Kit**



Item: RD01

Name: Dental Dam Kit

Packing: 1set/box,5boxes/ctn Carton Size: 450X270X305mm


### The Kit Contains:









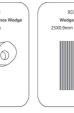



- Prevent instruments and materials from being inhaled, swallowed o damaging the mouth.
- Keep the operative field clean and dry, reach better performance of bonding materials.
- Isolate bacteria in saliva from splashing onto the tooth to avoid cros infection.
- Block lips and tongue for easy operation.
- Protect soft tissues and avoid injury from instruments and irritating
- Improved visibility of the operative site.
- Avoid the noise of the suction device.
- Save operation time, and improve treatment efficiency.

## **Punch**



- The punch head is combined with a spring and a tapered structure, which can be automatically adjusted to ensure the taper head is aligned with the center of the cutting edge when drilling.
- Hole diameter size: 0.8/1.0/1.4/1.6/1.8/2.0/2.3.

Item: RD04 Packing: 1 pc/box



## **Forceps**



- The forceps head with anti-shedding design, can prevent it from slipping when stretching the clamps.
- The right side of forceps is printed with positioning line. It can prevent extra stretching during operation, which will reduce resilience of clamps.
- Serrated design, accurate positioning, single hand operation, easy griping.

Item: RD05 Packing: 1 pc/box



#### **Frame**



- The frame is made of carbon fiber, with high hardness, strong toughness and is not easy to deform
- The shape is optimized to enhance the comfort when pulling the patient's
- The round edge can avoid skin damage, convenient to place the rubber dam.
- Can do heat sterilization.

Item: RD07/RD08

Packing: 1 pc/box,50boxes/ctn Remark: (RD07)L / (RD08)S



## **Wedge Floss**



- The wedge floss is with latex material and row design, the dentists can take appropriate quantity according to patients' gum gap.
- Excellent elasticity, in many cases, it can be used as a substitute for traditional clamps, which can quickly and easily fix the rubber dam.

Item: RD10

Packing: 1 set/box,68boxes/ctn Carton Size: 260X170X180mm



#### **Dental Dam**



- Malaysia imported latex, high tear tolerance, good flexibility.
- Blue: 6X6", Green: 5X5"

Item: RD02/RD03

Specification: (RD02)6X6" ,152X152mm (RD03)5X5",127X127mm

Packing: 1 pc/box,40boxes/ctn



## Clamp



Small size resin clamps for adult premolar large size for adult molar.

Item: RD11

Packing: (Z1+Z2)1 set/box,68boxes/ctn Carton Size: 260X170X180mm



The groove design of the clamp is convenient for the instrument to enter the groove and place the rubber dam. The two-wing design of the clamp fits firmly with the teeth. Using special material, the clamping force is stronger, and the matte surface design increases the texture of the product.

#### Clamp

















www.ajjhealthcare.com 25